

FOR HIGH RISK PATIENTS IN BOTH THE HOSPITAL AND COMMUNITY ENVIRONMENT


QUATTROTM overlay



The QUATTRO OVERLAYTM mattress system is designed for the prevention and treatment of pressure ulcers for high risk patients with limited mobility within both the hospital and community care environments.



The B.A.S.E.TM SEQUENTIAL cushion can be connected directly to the overlay for effective and affordable 24 hour pressure area care.

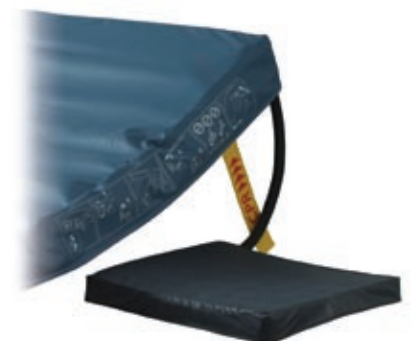

MICROPROCESSOR CONTROLLED PUMP OPERATES ON DEMAND ONLY




The QUATTRO OVERLAY™ features QUATTRO THERAPY™, a proven and effective method of preventing and treating pressure ulcers, unique to Talley Medical.

1 in 4 Cell Alternating Air Pressure Cycle provides pressure relief, lower average interface pressures and enhanced comfort.

ODS™ - Ortho-differential Support cell design provides optimum internal cell pressures for varying body shapes, as well as a feeling of security.

'Low Air Loss' Feature helps reduce moisture build-up, maintaining the patient's skin at an ambient temperature and lower humidity level.

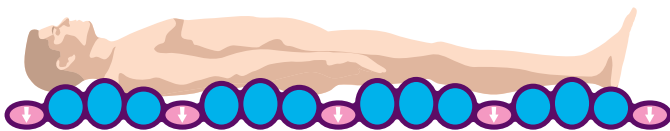

- PHASE FLOW™ pressure control pump operates on demand only, reducing power consumption and extending service life.
- Pump features adjustable patient comfort control, and visual and audible operation status and fault indicators.
- Simple and effective CPR facility.




QUATTRO OVERLAY™ MATTRESS



- Transport facility allowing unpowered support of the patient for up to 4 hours.

- Clinically proven for the prevention and treatment of tissue damage, the system provides effective therapy for patients with pressure ulcers up to Grade 3.

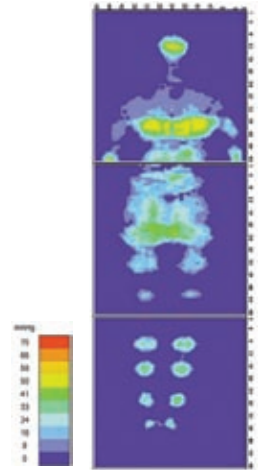

SPECIFICATION

Pump Dimensions	
Width	325mm
Height	233mm
Depth	165mm
Pump Weight	3.4kg
Cycle Type	1-in-4
Cycle Time	15 mins
Mattress Dimensions	
Length	1930mm
Width	855mm
Depth	130mm
Mattress Weight	6.3kg
Inner Construction	PU air cells
Cover Material	PU coated stretch nylon
Number of Cells	17
Max. User Weight Guide (kg/st)	160/25
Min. User Weight Guide (kg/st)	26/4

Talley Medical products are manufactured to comply with BSI, IEC, UL and other European safety standards. Talley Medical design and manufacture products to ISO9001, EN46001 and ISO13485 standards, and comply with Medical Directive (93/42/EEC). Every care has been taken to ensure that the information contained in this brochure was correct at the time of going to press. However, Talley Medical reserves the right to modify the specification of any product without prior notice in line with a policy of continual product development. Our standard terms and conditions apply. © Talley Group Limited 2003. All rights reserved.



TALLEY GROUP LIMITED
Premier Way, Abbey Park Industrial Estate,
Romsey, Hampshire SO51 9DQ England
Tel: +44 1794 503500 Fax: +44 1794 503555
e-mail: sales@talleygroup.com Web: www.talleygroup.com

