

The **FL-9000 Respiratory Humidifier** from Flexicare Medical meets the range of clinical needs in both hospital and home care environment with a single, multifunctional unit.


Simple to operate, with intuitive controls and easy to select functions, the FL-9000 Respiratory Humidifier offers a complete solution.

- Invasive or Non-invasive ventilation
- Heated wire or Non-Heated wire Applications



Multi-Function Display

Clear easy-to-read LED Display with the option to check Airway, Chamber or Heater Plate Temperature, or hold on a preferred parameter.

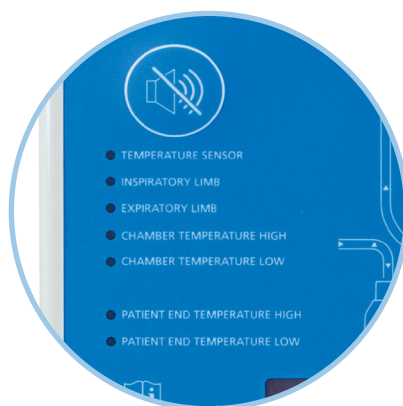

Alarm Indicator

Clear identification of fault source for quick clinical resolution.



Accessible

Easy-connect wire and probe sockets of different shape ensure correct connection.

Fully automatic for simple and rapid set-up, the FL-9000 Respiratory Humidifier has many advanced features, including preset default temperatures, a selectable temperature display, alarm panel with location indicators and a range of safety features.

Function Indicator

LEDs to indicate Function in Operation, Warning or Caution.

Warning System

Visual and audible alarms notify user of set-up errors, disconnection during use or delivered temperatures outside the set parameters.

Positive Connection

Heater Plate with spring loaded locking bar for secure location of humidification chamber.

Preventing Misconnection

Heated wire and temperature probe sockets designed to prevent misconnection.



Ordering Codes

FL9000	Respiratory Humidifier UK	230V with 3 Pin Plug	1
FL9000EU	Respiratory Humidifier European	230V with 2 Pin Plug	1
FL9000U	Respiratory Humidifier USA	110-120V with NEMA 5-15P 3 Pin Plug	1
FL9000AUS	Respiratory Humidifier Australia	230V with 3 Pin Plug	1
FL9000CH	Respiratory Humidifier China	220V with 3 Pin Plug	1
FL9000 INDIA	Respiratory Humidifier India	230V with 3 Pin Plug	1
FL9000-01	Single Heated Wire Humidifier Adapter		1
FL9000-02	Dual Heated Wire Humidifier Adapter		1
FL9000-03	Humidifier Temperature Probe (1.5m)		1



Accessories

FL2000-01	Rail Bracket Accessory	1
FL2000-02	C Clamp Accessory	1
FL2000-03	Mounting Bracket Accessory	1
FL2000-04	Mounting Kit Accessory	1
FL2000-05	Dual Mounting Kit Accessory	1

For optimal function and patient safety, the FL-9000 should only be used with breathing circuits, humidification chambers and other accessories that have been approved by Flexicare Medical.
Refer to the FL-9000 Respiratory Humidifier Operating Manual for further information on the technical specification, set-up and usage.